

NEW

CAV -47

VARIABLE FREQUENCY CAVITATION. DOUBLE CRYSTAL. SAFER and MORE EFFECTIVE.



LATEST TECHNOLOGY
PROFESSIONAL EQUIPMENT

- 2 BARIUM TITANATE CRYSTALS.
- VARIABLE FREQUENCY OF 35, 38 & 40 KHZ.
- PULSED and CONTINUOUS MODE.
- POWER CONTROL.
- FRIENDLY & EASY OF USE.
- SURGICAL STAINLESS STEEL HEADER.

CAV 47 - VARIABLE FREQUENCY CAVITATION

Cavitation is a novel body modeling technique that utilizes shock waves and pressure changes, forcing the adipocyte's implosion and its subsequent natural elimination.

CAV 47® with its double crystal system and shock wave power of 47 Watts grants higher application precision and stability.

CAV 47® - the best alternative to non-surgical liposuction - an outstanding option to fight localized adiposity assuring painless and short treatments with durable and visible results from the first session.

The therapist interpretation of the adipose panniculus will determine to set, via CAV 47®'s panel control, continuous or pulsed mode option.

With a variable frequency of 35, 38 & 40 KHz, CAV 47® grants a friendly and easy way of managing any adipose tissue deepness and consistency.

WELCOME TO TECHNOLOGY! WELCOME TO NOTO!

NOTO
THE DIFFERENCE OF BEING UNIQUE

NOTO proposes a different course for aesthetics, a new scenario where therapies combination, cutting-edge technology, and quality products and supplies are fundamental for visible, effective, and durable results.

THE BEST ALTERNATIVE

Design for body modeling and non-invasive, safe, and effective fat elimination. CAV 47® is the best option to fight localized adiposity and collaborate with cellulite improvement, assuring visible and durable results from the first sessions.

CAV 47® creates the implosion of micro bubbles caused by pressure changes disintegrating localized fat cells. These turn into a liquid substance which is removed through the lymphatic system and urinary tract.

TREATMENTS.

- Pre & post liposuction.
- Modeling of body outline.
- Hard areas to slim down.
- Improve cellulite external aspect.


AREAS TO TREAT

- Inner and external leg sides
- Arms
- Abdomen (left, right, and central regions)
- Glutei.
- Orange peel cellulite.

CAV-47 IPX 0 DEVICE

MANUFACTURED BY:	NOTO GROUP S.A.
COUNTRY:	ARGENTINA
NAME OF PRODUCT:	CAV-47
WIDE:	405 mm.
LENGTH:	310 MM.
HEIGHT:	160 MM.
WEIGHT:	4.6 KG.
POWER SUPPLY VOLTAGE:	98 to 240 VAC 50/60HZ
ENERGY CONSUMPTION:	45V
POWER:	47 WATT
PROGRAMMING:	MICROCONTROLLER
FREQUENCY:	35, 38 & 40 KHZ +/- 10%
DISPLAY:	128X64 GRAPHIC
SOURCE:	INTERNAL SWITCHING
TECHNOLOGY:	SMT
TIMER:	MANUAL
ENTRANCE FUSE TYPE:	F 1.5 AL 250 V
EMITTER ENTRANCE FUSE:	F 1 AL 250 V

NOTE: DEVICE SPECIFICATIONS MAY CHANGE WITHOUT PREVIOUS NOTICE.

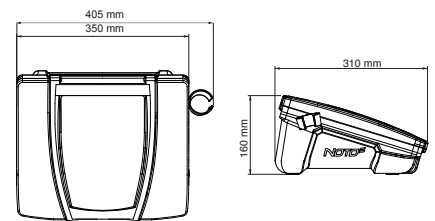

CAV-47

NOTO supports all clients with a set of tools available at all times.

Tools as training, post-sale service support, and assured spare parts. All these instruments with the aim of making Noto's products your best choice.

More equipments by NOTO:

TRIMAX - Resistive Radiofrequency of 3 emitters and cross fields.

N 1104 -Biphasic Electrotherapy. T.E.N.S. + RUSSIAN WAVES+ SQUARE WAVES. 26 Programs.